


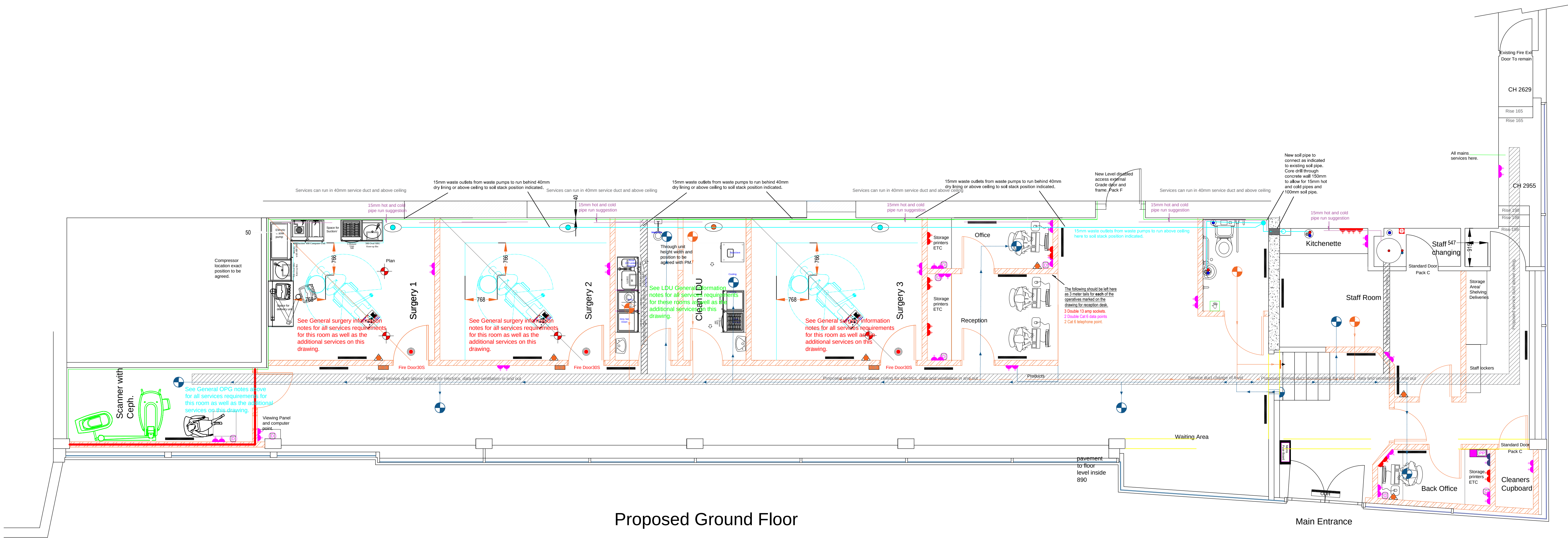

Proposed Ground Floor

IMPORTANT NOTES TO CONTRACTOR:

♦ **SUCTION & WASTE:**
All waste and suction pipe work to be no less than 40mm pipe.
Suction pipe work must be UPVC/Polypropylene material. **DO NOT USE ABS PLASTIC PIPE.** Sweep bends are only to be used and the number of bends must be kept to a minimum (preferably no more than 2 or 3 in any run). All waste must be gravity wastes.
WASTES MUST HAVE A TRAP FITTED IN LINE TO PREVENT UNWANTED ODOURS TRAVELING BACK TO SURGERY. TRAP MUST BE EASILY ACCESSIBLE.
We also request that the contractor installing the pipe work indicates on the actual pipe work which pipe is suction, waste & vent pipe (if applicable).

♦ **WRAS' approved Double Check Valve** should be fitted in a suitable serviceable position in the water line supplying the dental treatment centre.

♦ **POWER:**
When a power supply is less than 1.5mtrs away from a water supply an R.C.D. (residual current device) must be fitted.

♦ **ISOLATION:**
The mechanical services contractor must ensure that the water and air supplies to each individual surgery can be isolated to enable the dentist to turn off the water & air supplies at the end of surgery hours. These isolation valves must be fitted in an easily accessible position. It is also the responsibility of the mechanical services contractor to supply & fit the relevant parts for this isolation.

♦ **COMPRESSED AIR:**
It is recommended that air lines are run in 15mm copper piping. An air/water separator should be fitted in line at the pipe entry point into the surgery. The pipe should then be lagged from the air/water separator to the pipe termination point. Adequate earth bonding will be required to the floor box.

♦ **FINAL CONNECTIONS TO CABINETRY:**
Cabinetry fitters will not make the final connections to the cabinetry services (sink services & electrics). These final connections must be made by the contractors.

WHEN ROUTING THE REQUIRED PIPE WORK, PLEASE BE AWARE THAT SUBSTANTIAL FIXING BOLTS MAY BE USED WHEN SECURING THE CHAIR TO THE FLOOR. PLEASE CHECK WITH THE SPECIALIST FOR THE APPROXIMATE POSITION OF THE FIXING POINTS. SO THAT THE PIPE WORK CAN BE RUN TO AVOID THESE POSITIONS.

General Dental Chair Service Requirements For Each Surgery:

Exact chair position to be confirmed
Either at the chair box or within the cabinets please allow for mains cold water, compressed air, electric, suction pipe, waste & electric switch wires from suction motor.

General Service Requirements For Each Surgery:

Six double sockets, three data points, power for suction motor, mains cold water, hot water and wastes for two sinks as well as that shown on the drawing.

General OPG X-Ray Requirements

Wiring (power) for x-ray and wiring for the controls
Wiring for the wall mounted timer / operating switch.
A wall mounted switched fused spur to isolate & power the x-ray.

A PC will be necessary for the OPG linked to the machine and patch panel. Provide all power and data services for this. Exact PC position to be agreed. Specialist to forward detailed layout of room.

Client Information

Any indications given in respect of radiation protection, the siting of x-ray equipment, and / or x-ray controls is for guidance purposes only. It is for your appointed equipment specialist and RPA to review the environment in which the x-ray machine is to operate and to ensure your compliance under IRMER. This drawing should not be used as the basis for any design, construction work or the provisioning of services involving x-ray equipment without the approval of your RPA and chosen equipment specialist.

RPA

Any information that Archintelligence has provided regarding the type of X-Ray wall protection, has been given to assist the RPA in their design report for the client and the contractor for cost purposes. This information has been given from our best judgement and at no stage has it been reviewed, or supported by an RPA specialist prior to contract tendering. Therefore under no circumstances will Archintelligence accept responsibility for any inaccuracies.

Clean LDU

Allow to

Ventilation
See service drawing

Mechanical ventilation provided to WC capable of extracting 10 l/s per person.
Allow for fresh in via vents shown and extract out to the LDU as per the service drawing. All to be ducted to outside, air and flow rate commissioning certificates to be provided on completion.

NOTE: There is limited access to the external due to shop front glazing and low ceilings on car park outside. Suggested route is above existing flat roof. Flat roof itself could be used or each room has drilled access through to car park side (if landlord allows-PM to advise).

Air admittance valves Durgoo or similar to WC should terminate above the highest waste appliance.

Lighting and Power General

New lighting to achieve a luminous efficiency of 60 lamp lumens per circuit watt. 75% of all new lighting to be LED's. New control sockets to be 450-1200mm above floor level. See lighting drawing.

Electrics

Electrics to be installed by a competent person in all accordance with BS7671:2008 Part P. Commissioning certificate and test commission certificates to be provided no later than 30 days after completion of work.

Hot Water and Heating

The existing hot water system is to be retained to supply all hot water to the practice. See service drawing for outlets.
Existing heating system is also to be adjusted to accommodate the practice layout. Correct BTU's to room space to be provided.

Clean And Dry LDU General Requirements
Allow for hot, cold and waste for a total of 4 sinks, 10 double sockets, 2 data points 2 phone points.


Services Legend- Exact Service Positions TBA

Symbols	Description	Symbols	Description
	Proposed Double Socket 350 H		Overhead Door Heater
	Proposed Double Socket 875 H		Ceiling mounted mechanical extract fan-to be ducted to external walls as per the drawing. Must have louvred external cover with fly-screen.
	Proposed Double Socket for patch panel		Suggested mechanical extract pipe runs. Must have louvred external cover with fly-screen.
	Proposed Single 13 AMP Switched fused spur wall mounted 2290 H for extract.		Electric energy efficient 450 W (with Eco benefits such as Haverland). Slim line radiators, exact locations to be agreed on site. BTU to room space must be correct.
	Proposed Single 32 AMP Switched fused spur wall mounted 2290 H for air Megaflo.		Wall mounted grid for fresh air In-to be ducted through external wall. Does not have to be mechanical. Must have louvred external cover with fly-screen.
	Switch Bank Isolator for all power to surgery except CPU unit and fridge if installed.		Ceiling mounted grid for fresh air In-to be ducted through external wall. Does not have to be mechanical. Must have louvred external cover with fly-screen.
	Proposed Phone outlet exact position to be agreed.		Suggested fresh air pipe runs. Must have louvred external cover with fly-screen.
	Proposed Wall Mounted Hand Dryer		Suggested extract air pipe runs. Must have louvred external cover with fly-screen.
	Proposed Single Cat 6 Data Socket wall mounted at 1200mm high to centre of socket.		Waste water pump Saniflo or similar connected to suction waste and sinks waste. Exact position to be agreed with PM. Allow 240 power supply for same.
	Proposed Single Cat 6 Data Socket wall mounted to below work top height (350) unless otherwise agreed		Waste water pump Saniflo or similar connected to sinks waste. Exact position to be agreed with PM. Allow 240 power supply for same.
	Proposed Single Cat 6 Data Socket ceiling mounted.		
	Proposed Patch Panel exact position to be agreed.		
	Proposed Single Cat 6 Data Socket at 750mm to centre of socket off finished floor.		
	Fuse board at high level exact position to be agreed.		
	Proposed cold water for WC, drinks machine or water tank.		
	Proposed hot/cold and waste exact position to be agreed.		
	1 x Megaflo to serve all building with hot water.		

Do not scale off drawings always use site measurements.

Revision Register			
No.	Description	Date	Approved by
2	Compressor position changed	27-8-2019	

Note:

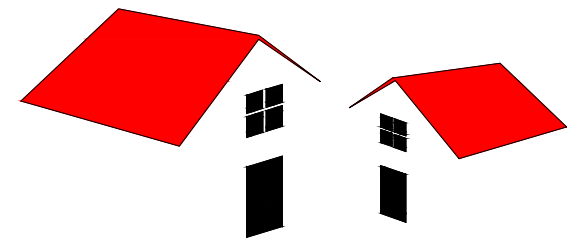

Archintelligence

www.archintelligence.com

E: info@archintelligence.com

Tennysen House
Cambridge Business Park
Cambridge CB4 0WZ

CLIENT:	Set and Roopa Regarding The New Practice at Lammore House 370 High Road Wembley London HA9 6AX	DATE:	2-9-2019	Approved by:	DATE:
DWG NO:	SARLH 2-9-19 PSD	SCALE:	1:50 at A1		
DWG TITLE:	Proposed Service Drawing	DRAWN BY:	RG		